



Full automation  
through  
consolidation

Automated Urine Sediment Analyzer

AUTION EYE  
AI-4510

Fully Automated Urine Analyzer

AUTION MAX  
AX-4060



Automated Urine Sediment Analyzer

# AUTION EYE

## AI-4510

Automatic Dilution Function  
Atlas Image Collection  
microscopic imageautomatic



Compact Footprint  
Understanding the value shift  
Gating FunctionDilution Function

Fully Automated Urine Analyzer

# AUTION MAX

## AX-4060

## Fewer Manual Steps

### Automatic Dilution Function

Autodilution of highly concentrated samples  
Reduces the need for further testing

## Unique Training Tool

### Atlas Image Collection

12-category library  
New, interesting images can be added to the collection

## Sensory identification

### Onscreen image is presented as a (synthetic) microscopic image

Colour images assist in particle identification  
Summary screen allows rapid assessment of particle frequency and size

## Ease of implementation

### Compact Footprint

AUTION EYE can be connected to the AUTION MAX series (Smart work cell)

## High Quality Results


### Understanding the value shift

Compare chemistry test results with sediment analysis (Cross-Check Function)  
Improves reliability of results

## Reduces need for further testing

### Gating Function


Urinalysis results of AX will determine if further sediment analysis is necessary  
User can define testing rules

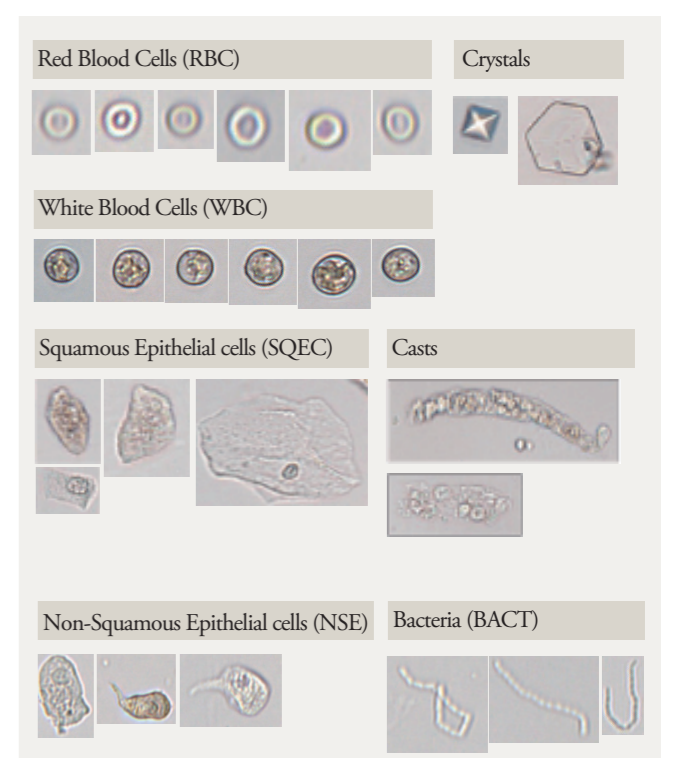

Example of particle images

# Atlas Image Collection (reference figures)

Atlas includes preinstalled images. Users can add their own images from measurement results. Atlas is a useful reference for help with manual classification.

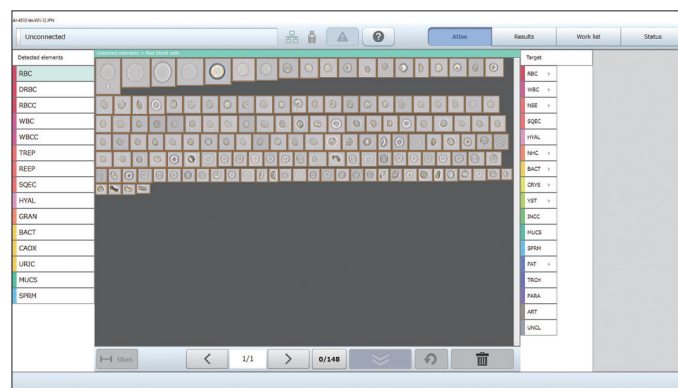


## Summary Screen

Formed element images are displayed as a (synthetic) microscopic image.

The formed element selected in the list is highlighted.

Click on the mouse wheel for image selection.

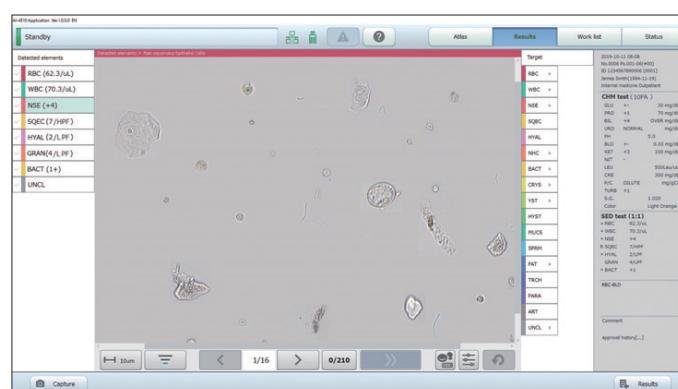

## Related products



New reagent

### Uriflet S 11UA

Uriflet S 11UA strips can help with the screening and diagnosis of diabetic nephropathy and CKD (chronic kidney disease).

Asymptomatic early stage renal disease can be identified accurately by simultaneous microalbumin/creatinine measurement, allowing for greatly improved prognosis.

Test strip	Measurement items											Calculated items	
	GLU	PRO	BIL	URO	pH	BLD	KET	NIT	LEU	CRE	ALB	P/C	A/C
11UA	●	●	●	●	●	●	●	●	●	●	●	●	●



Semi-automatic Urine Analyzer

### Aution|Daten

AE-4070

Compact design:

206(W)×365(D)×180(H)mm/5.0 kg

High throughput:

514 tests/hour (7 seconds/test)

Auto-start function:

Measurement starts automatically when a test strip is set onto the tray.

A touchscreen allows for easy and intuitive operation The successor of Aution Eleven.

This can be backup for AX-4060.

## Legal manufacturer

**arkray factory, inc.**

1480 Koji Konan-cho, Koka, Shiga-ken 520-3306, JAPAN

## European representative

**arkray europe, b.v.**

Prof. J.H. Bavincklaan 2, 1183 AT Amstelveen, THE NETHERLANDS

<https://www.arkray.eu/english/>

\*Designs and specifications may be changed without prior notice.

**arkray global business, inc.**

Yousuien-nai, 59 Gansuin-cho, Kamigyo-ku,

Kyoto 602-0008, JAPAN

TEL +81-75-662-8979 FAX +81-75-431-1202